



Rotary Evaporator Chiller




Unique Features

- Digital display for live working status of the system.
- Inbuilt smart technology with emergency alert and auto cut off for low water flow and low water level.
- User friendly interface for temperature settings.
- RS 232 communication with IOT and data login.
- We specialize in making customized chillers with heat dissipation capacity from 1 KW to 20 KW.







Home Screen

Home screen or front page on the display consists of, Flow rate in L/min, Water line pressure in kg/cm and Tank water level in percentage(%) .

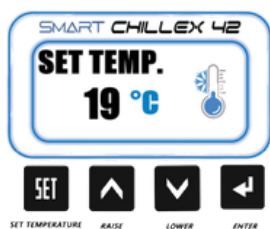



Operation & Buttons

On Pressing set button. The display moves to second screen where you can set temperature using the RAISE button to increase setpoint, LOWER to decrease the setpoint and ENTER button to send value to the system.

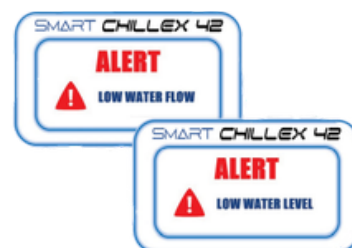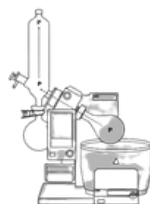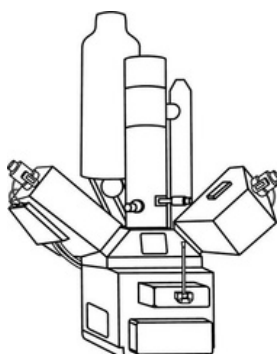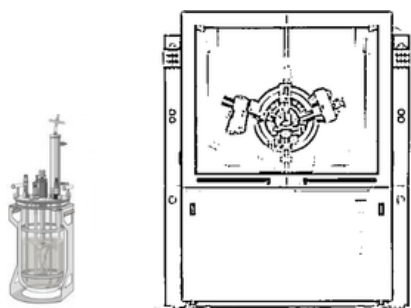

Automated Alerts

Error screen automatically pops up when system has Lower Water Flow, Low Water Level, saving the system from dry run damage.

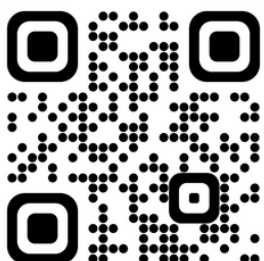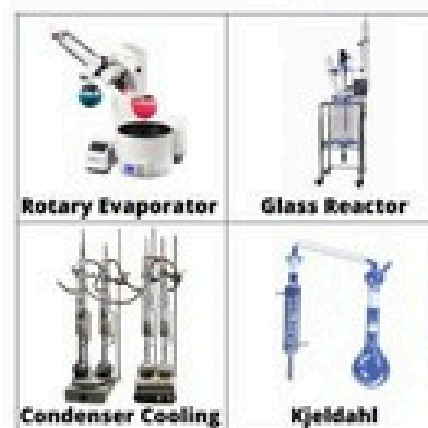

Hexon Instruments Pvt Ltd

EL-91, Electronic Zone, MIDC Industrial Area,
Mahape, Navi Mumbai, Maharashtra 400710



+918657448273



info@hexoninstruments.com
sales@hexoninstruments.com